

Instruction Manual

Axillary Crutch Grey



CTALGRY

SectorCare

<https://sectorcare.com.au>



Product Installation



01 Handle

02 Grip


03 Height Adjustment

04 Anti-Slip Tip

Parts List:



Metal Underarm Tube*2

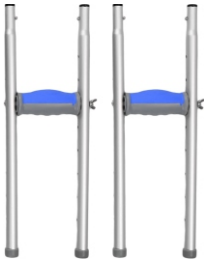

H-shaped Round Tube*2



Grey Adjustable Tube*2



Metal Leg Tube*2

1. Take out all parts and accessories after opening the package.



2. Align the connecting joints of the underarm tube and H-shaped round tube, insert into place and secure.



3. Align the connecting joints of the leg tube and grey adjustable tube, insert into place and secure.



4. Finally, align the connecting joints of the two pre-assembled components, insert into place and secure.



LEGAL STATEMENT

This manual and its contents are protected under the Copyright Act 1968 (Cth). No part of this publication may be reproduced, distributed, modified, stored in a retrieval system, or transmitted in any form or by any means without the prior written permission of FRG Pty Ltd. Unauthorized use is strictly prohibited.