

Instruction Manual

Shower Chair with Armrest and 360 degrees



SBST36001MSIL

SectorCare

<https://sectorcare.com.au>



Product Overview



1 Back Panel

2 Hand Grip

3 Seat Board

4 Metal Base Support

5 Foot Tube

6 Anti-slip
Rubber Tip

Product Installation

Parts List:



Seat Board*1



Back Panel*1



Hand Grip*2



Foot Tube*4



H-shaped Round Tube*2




Clip*1

Product Installation


Take out all parts and accessories from the carton.



Take out the H-shaped round tubes, align them with the connectors on both sides of the seat main frame respectively, insert and fix, as shown in the figure.

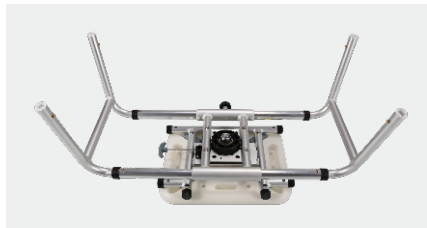

Take out four leg tubes, insert them into the four corresponding connectors of the main frame respectively and secure, as shown in the figure.

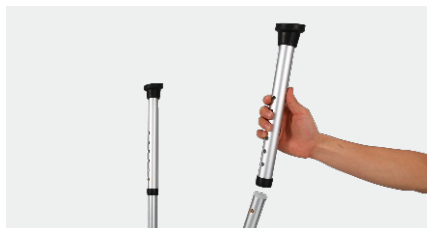

Align the armrests with the connectors on both sides under the seat board respectively, insert and secure them, as shown in the figure.



Take out the backrest, place it at the connector under the seat board, insert and fix it, as shown in the figure.



The buckles can be installed on the leg tubes as needed (the buckles are for holding the shower head), as shown in the figure.



LEGAL STATEMENT

This manual and its contents are protected under the Copyright Act 1968 (Cth). No part of this publication may be reproduced, distributed, modified, stored in a retrieval system, or transmitted in any form or by any means without the prior written permission of FRG Pty Ltd. Unauthorized use is strictly prohibited.

© 2026 FRG Pty Ltd. All rights reserved.